


T4 LINE / Air Motor

sirona
The Dental Company

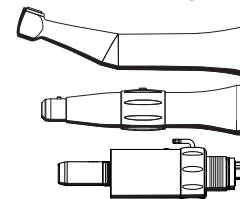

English

Operating Instructions

Table of Contents

1	Before you begin ...	5
1.1	Structure of the document	7
1.2	Service life of Sirona instruments	10
2	Safety information	11
3	Technical description	13
3.1	Task	13
3.2	Assembly of the T4 LINE	14
3.3	Design of Air Motor	15
3.4	Technical data	16
4	Preparation	20
4.1	Initial start-up and longer breaks in use	20

Operating Instructions T4 LINE / Air Motor


Table of Contents

4.2	Prior to starting the work day.....	20
4.3	Before each patient	20
5	Operation	21
5.1	Connect motor to the supply hose.....	22
5.2	Replacing the instrument.....	23
5.3	Inserting and removing burs.....	23
5.4	Inserting and removing the disposable prophylaxis head	27
5.5	Connect and remove external spray	29
5.6	Setting drill direction on the motor.....	30
6	Follow-up	31
6.1	After each treatment session.....	31
6.2	At the end of the work day.....	32
6.3	Once per work week.....	32


Table of Contents

Operating Instructions T4 LINE / Air Motor



7	Conditioning	33
7.1	Conduct pre-disinfection.....	33
7.2	Automated cleaning and disinfecting	34
7.3	Manual cleaning and disinfection	37
7.4	Manual maintenance	39
7.5	Sterilizing.....	43
8	Maintenance	46
8.1	Replacing the sealing washer	46
8.2	Lubricate the lock washer.....	47
9	Spare parts and consumables	48
10	Storage and transport conditions	49
11	Disposal	50

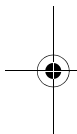

Intended use

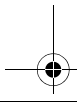

63 97 959 D3451.201.06.06.02

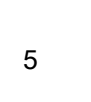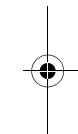

1 Before you begin ...


T4 Line / Air Motor complies with state-of-the-art technical regulations.
T4 Line / Air Motor complies with the ISO 14457 standard.

The motor is to be connected to an interface according to ISO 9168.

1. Read the operating instructions prior to using T4 Line / Air Motor.
2. Only use T4 Line / Air Motor for the applications described in the operating instructions.
3. Observe the applicable hygiene standards, occupational safety regulations and accident prevention measures for T4 Line / Air Motor.

The instrument is used to hold and drive burs for the purposes of rotary processing and is intended for the following uses in dentistry:

- Removal of carious material and fillings; cavity and crown preparation of hard tooth substance; processing tooth and

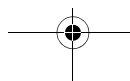




6397959.book Seite 6 Mittwoch, 26. Juli 2017 2:52 14




1 Before you begin ...

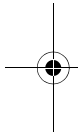

Contraindications
Target group



6

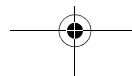


Operating Instructions T4 LINE / Air Motor

- restoration surfaces; processing dental prostheses: crowns, fillings, bridges (contra-angle handpiece)
- Prophylaxis (contra-angle handpiece)
 - Endodontics (contra-angle handpiece)
 - Implantology, surgery, processing bones (contra-angle handpiece)
 - Prosthetics (straight and contra-angle handpiece)


The micro-motor serves as a drive for rotating and oscillating straight and contra-angle handpieces in connection with a dental treatment center.

None

This product is intended only for use by trained dental personnel in dental practices and laboratories.

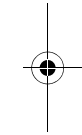


63 97 959 D3451.201.06.06.02

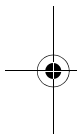

Warnings

Instructions for use

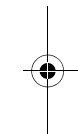

1.1 Structure of the document

1.1.1 Labeling of information

➤ To prevent injuries, please observe the warnings.

Warnings are labeled as follows:

⚠ **DANGER!** indicates a danger **leading** to death or serious injury if not avoided.

⚠ **WARNING!** indicates a danger that **may lead** to death or serious injury if not avoided.

⚠ **CAUTION!** indicates a danger that **may lead** to injury if not avoided.

➤ To prevent material damage and additional expenses, please observe all instructions for use.

Instructions for use are labeled as follows:


NOTICE! indicates measures for the prevention of material damage.




6397959.book Seite 8 Mittwoch, 26. Juli 2017 2:52 14



1 Before you begin ...

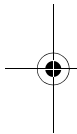


8

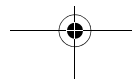

Operating Instructions T4 LINE / Air Motor

IMPORTANT: indicates information on the avoidance of additional expenses and other important information.

Tip: indicates information for facilitating work.

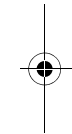

63 97 959 D3451.201.06.06.02




1.1.2 Formats and symbols

The formats and symbols used in this document have the following meaning:

✓ Prerequisite 1. First action step 2. Second action step or > Alternative action ↪ Result ➤ Individual action step	Requests you to do something.
Use of formats and symbols [→ 9].	Identifies a reference to another text passage and specifies its page number.
• List	Identifies a list.




1 Before you begin ...

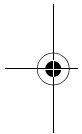

10



Operating Instructions T4 LINE / Air Motor

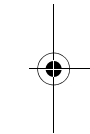

1.1.3 Abbreviations/codes used

CA	Contra-angle handpiece bur instrument
HP	Straight handpiece bur instrument

1.2 Service life of Sirona instruments

When used as intended:

- Non-moving parts have a typical service life of 3 years
- Moving parts have a typical service life of 1.5 years

No warranty claim can be inferred here, as wear may occur earlier or later than indicated above depending on use, frequency of sterilization, and frequency of maintenance.

Elastomers, e.g. O-rings, must be replaced depending on their degree of wear.

For safety and technical reasons, check the clamping system of the contra-angle and straight handpiece burs on an annual basis.

63 97 959 D3451.201.06.06.02

Obligations of the user

Preventing the spread of infections and cross contamination

Instrument head overheating

Hot motor

63 97 959 D3451.201.06.06.02

2 Safety information

- Use only fault-free materials that do **not** deviate from the specified data [→ 16].
- Protect yourself, patients, and others against any foreseeable dangers. To do this, follow the safety information.
- Comply with the Intended use of the equipment.
- You should always keep these operating instructions within reach for further reference.

Prevent the spread of infections and cross contamination between patients, users, and third parties. Sterilize equipment after each patient.

Take the appropriate hygiene measures, e.g. wear protective gloves.

If the instrument is defective, the area around the instrument head may heat up, thus creating a risk of burning the patient's oral mucosa.

If the motor is defect, it may get hot.



2 Safety information

Malfunction or damage

Potential explosion hazard

Repair

Spare and accessory parts

Operating conditions:



Operating Instructions T4 LINE / Air Motor

Discontinue use immediately in case of malfunction, unusual or different sounds or damage. Damaged instruments may cause injury. Notify the dental dealer or the manufacturer.

Do **not** use this product in areas subject to explosion hazards.



Do **not** repair the instrument yourself.

Use only original Sirona parts or parts approved by Sirona. Safe operation is **not** guaranteed for parts that have not been approved by Sirona.

Temperature: +10 °C – +30 °C

If you have any questions, please contact your dental depot or the manufacturer.



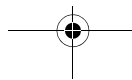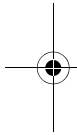


3 Technical description

3.1 Task

The instrument transmits the driving power and speed of the electric motor or air motor to the preparation tool.

The motor converts air to rotation and transfers rotation to the instrument.

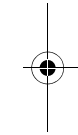


3 Technical description

Operating Instructions T4 LINE / Air Motor


3.2 Assembly of the T4 LINE

Contra-angle handpiece

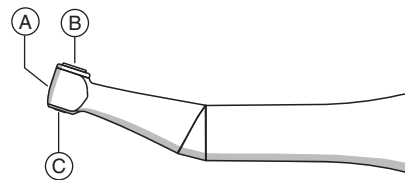


A	Instrument head
B	Push button
C	Opening for the clamping system

Straight handpiece

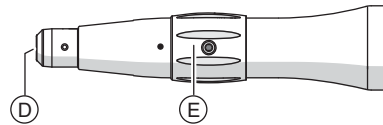

D	Opening for the clamping system
E	Center part of handpiece

Operating Instructions T4 LINE / Air Motor

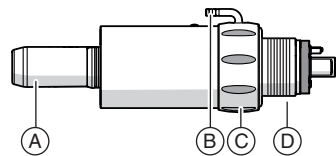

3.3 Design of Air Motor

A	Handpiece holder
B	External spray connector
C	Drill direction switching ring
D	Coupling connection (here: Midwest/ISO 4-hole)

3 Technical description

3.4 Technical data

Instrument

	T4 LINE B 40	T4 LINE BH 40
Gear ratio	~ 1:1	~ 1:1
Maximum operating speed in rpm	~ 40000	~ 40000
Clamping system	CA	HP/CA
Instrument coupling	INTRAmatic®	INTRAmatic®

Air motor

	Air motor M / Air motor B
Motor speed in rpm	> 18 000
Torque in N/cm	~ 2.1
External spray function	x
Handpiece holder – ISO 3964	x

Bur

	CA	HP
Shank diameter in mm	2.334 - 2.350	2.334 - 2.350
Maximum total length in mm	34	50

3 Technical description

Operating Instructions T4 LINE / Air Motor

	CA	HP
Maximum working diameter in mm	-	-
Gripping length in mm	≥ 11	≥ 15
Standard	ISO 1797-1, Type 1	ISO 1797-1, Type 2

Dental treatment center

	T4 Line / Air Motor
Driving air pressure in bar (flowing, dry, and clean)	2.3
Return air pressure in bar	< 0.15
Air consumption in NI/min	~ 50

Operating Instructions T4 LINE / Air Motor



63 97 959 D3451.201.06.06.02

3 Technical description

	T4 Line / Air Motor
Spray air pressure in bar	2.3 ± 0.2
Spray water pressure in bar	2.0 ± 0.2
Maximum water temperature in °C	40
Recommended water content of spray	> 50 ml/min

This product bears the CE mark in accordance with the provisions of the Council Directive 93/42/EEC of June 14, 1993 concerning medical devices (MDD).

4 Preparation

Operating Instructions T4 LINE / Air Motor

4 Preparation

4.1 Initial start-up and longer breaks in use

- Sterilize T4 Line / Air Motor and accessories prior to startup.
- Clean and maintain T4 Line / Air Motor after longer breaks in use.

4.2 Prior to starting the work day

- Purge the water and air channels for 30 seconds.

4.3 Before each patient



⚠ CAUTION! Insufficient cooling leads to overheating of the preparation site and damage to the tooth substance. Ensure that the water content is > 50ml/min.



6397959.book Seite 21 Mittwoch, 26. Juli 2017 2:52 14



Operating Instructions T4 LINE / Air Motor

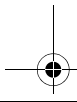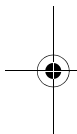

63 97 959 D3451.201.06.06.02



5 Operation


NOTICE! Use only burs and diamond polishers that are sharp and undamaged. Use clean burs and diamond polishers to avoid dirt in the clamping system.

⚠ CAUTION! The handpiece must only be operated with the clamping system closed.



⚠ CAUTION! Insufficient cooling leads to overheating of the preparation site and damage to the tooth substance. Ensure that the water content is > 50 ml/min.

⚠ CAUTION! Do not pull the patient's cheek back with the contra-angle handpiece while the motor is running. This would actuate the pushbutton, thus creating a risk of burning the patient's oral mucosa.

IMPORTANT: For application details and operating data of the different preparation tools, please refer to the information provided by the manufacturer.

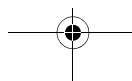


5 Operation

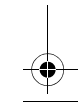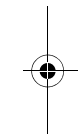

21



5 Operation

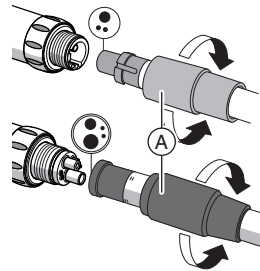

22

Operating Instructions T4 LINE / Air Motor

IMPORTANT: Let the handpiece idle at least 30 seconds before use.

5.1 Connect motor to the supply hose

1. Attach the motor to the supply hose, noting the various tube diameters.
2. Screw the screw cap (A) hand-tight. Hold the motor securely when doing this.

Does water leak out between the motor and the hose connection?

1. Detach the motor from the supply hose.
2. Reconnect the motor to the supply hose. Make sure it is connected properly.
3. If water still leaks, replace the sealing washer [→ 46].

63 97 959 D3451.201.06.06.02

5.2 Replacing the instrument

⚠ CAUTION! The instrument should only be fitted or removed when the motor is at standstill.

NOTICE! Do not use instruments with fiber optics on motors that do not have fiber optics.

Attaching the instrument

- ✓ The motor is at a standstill.
- Insert the instrument until it snaps into place.


Removing the instrument

- ✓ The motor is at a standstill.
- Detach the instrument. Do not pull on the supply hose while doing this.

5.3 Inserting and removing burs

IMPORTANT: Check the push button to make sure it moves freely!

5 Operation

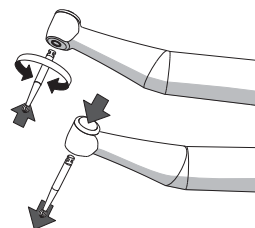

24

Operating Instructions T4 LINE / Air Motor

NOTICE! Use only burs and diamond polishers that are sharp and undamaged. Use clean burs and diamond polishers to avoid dirt in the clamping system.

⚠ CAUTION! The handpiece must only be operated with the clamping system closed.

Inserting the contra-angle handpiece bur

✓ The motor is at a standstill.

1. Insert the bur without pushing the button.
2. Turn the bur back and forth gently until it snaps into place.
3. Pull and turn the bur to check that it is firmly in place.


Removing the contra-angle handpiece bur

✓ The bur must not be moving.

➤ Press the button and remove the bur.

63 97 959 D3451.201.06.06.02

Operating Instructions T4 LINE / Air Motor

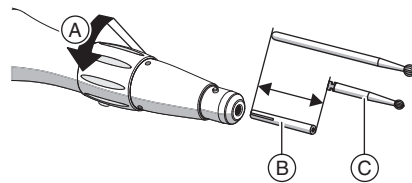


Inserting the straight handpiece bur

✓ The motor is at a standstill.

1. Turn the center part of the handpiece in the direction of the arrow (A) until it reaches the stop.
2. If you are using a **contra-angle handpiece burr instrument (C)**: insert the pin (B) into the handpiece, slit end first. The pin compensates for the difference in length between straight and contra-angle handpiece burr instruments.
3. Insert the burr instrument until it reaches the stop.

5 Operation


5 Operation

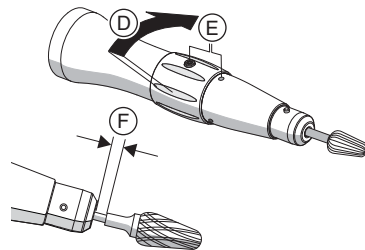

Operating Instructions T4 LINE / Air Motor

4. Clamp the burr instrument by turning the center part of the handpiece in direction (D) until it reaches the stop.
 - ↳ When the screw and the marks (E) line up with one another, the burr instrument is clamped.
 - ↳ When using handpiece burr instruments with a reduced shank: always make sure that the reduced shank does not come into contact with the opening for the chuck (F).
5. Pull and turn the burr instrument to check that it is firmly seated.

Operating Instructions T4 LINE / Air Motor

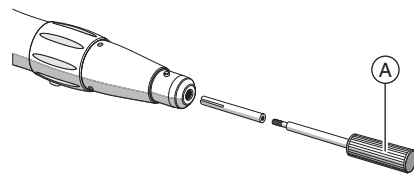

Removing the mandrel adapter

1. Screw in the supplied tool (A).
2. Open the clamping system.
3. Pull out the mandrel adapter.


5.4 Inserting and removing the disposable prophylaxis head

The handpiece fits disposable prophylaxis heads (e.g., Doriot type made of plastic).

⚠ CAUTION! The handpiece may be operated only with the clamping system closed.

IMPORTANT: Reduce the maximum operating speed of the handpiece according to the manufacturer's specifications on the dental unit

5 Operation

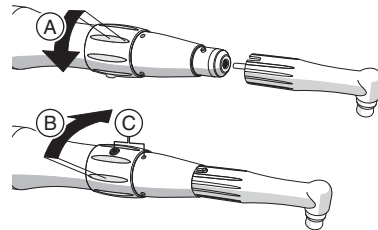

28

Operating Instructions T4 LINE / Air Motor

Inserting a disposable prophylaxis head

✓ The motor is at a standstill.

1. Turn the center part of the handpiece in the direction of the arrow (A) until it reaches the stop.
2. Insert the drive shaft of the disposable prophylaxis head into the handpiece until it reaches the stop.
3. Close the clamping system by turning the center part of the handpiece in the direction of the arrow (B) until it reaches the stop.


✎ When the screw and the marking (C) are lined up with one another, the disposable prophylaxis head is clamped.

Removing a disposable prophylaxis head

✓ The motor is at a standstill.

1. Open the clamping system.
2. Remove the disposable prophylaxis head.

63 97 959 D3451.201.06.06.02

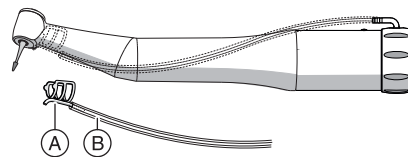

5.5 Connect and remove external spray

Connect external spray

1. Attach the instrument to the motor.
2. Attach the silicone hose (B) to the small tube on the spray clip (A) and motor.
3. Attach the spray clip to the instrument.

Remove external spray

1. Detach the silicone hose from the motor.
2. Remove the spray clip from the instrument.
3. Flush the spray clip and silicone hose with clear water.

5.6 Setting drill direction on the motor

The motor has two drill directions:

- Clockwise (forward), marked F
- Counterclockwise (reverse), marked R


Set forward




➤ Turn the switching ring clockwise to the stop.

Set reverse

➤ Turn the switching ring counterclockwise to the stop.

IMPORTANT! If the external spray connection is in the middle between R and F, the motor is switched off.

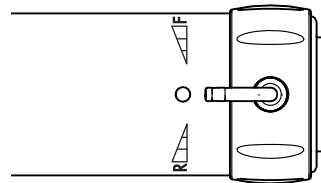

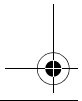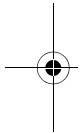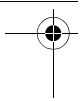

6 Follow-up

6.1 After each treatment session

NOTICE! Condition immediately, or at the latest, one hour after treatment.





NOTICE! Only use Sirona T1 Spray.

NOTICE! The air motor is **not** suitable for use with the Sirona DAC UNIVERSAL or with a cleaning and disinfection device.

✓ The bur must not be moving.

✓ Wear appropriate protective clothing.

1. Purge the water and air channels on the treatment center for 30 seconds.
2. Remove the bur with tweezers.
3. Pre-disinfect directly at the treatment center [→ 33].
4. Remove T4 Line / Air Motor.

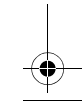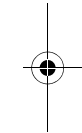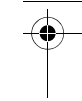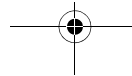

6 Follow-up

Operating Instructions T4 LINE / Air Motor



5. Transport T4 Line / Air Motor to the hygiene room in a suitable transport container.
6. Conduct automatic conditioning [→ 34]. Manual conditioning [→ 37] is possible in exceptional cases if the national/local regulations are followed.
7. Apply spray [→ 39] to T4 Line / Air Motor.
8. Sterilize T4 Line / Air Motor and the accessories [→ 43].

6.2 At the end of the work day

- Apply spray to the instrument [→ 39].

6.3 Once per work week

- ✓ The motor has been conditioned.
- Lubricate the lock washer [→ 47].

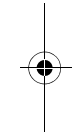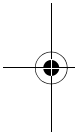

7 Conditioning


7.1 Conduct pre-disinfection

- ✓ Wear appropriate protective clothing.
- ✓ All disinfectants must be approved in your country and have proven bactericidal, fungicidal and virucidal properties. Only use disinfectants with **no** protein-fixing properties.

1. Spray the surface with disinfectant.
2. Wipe the disinfectant away with a cloth.

👉 For further processing, T4 Line / Air Motor should be dry and free of residue.

Please observe the manufacturer's instructions for using instrument disinfectants.

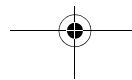

7.2 Automated cleaning and disinfecting ...


7.2.1 ... with a DAC UNIVERSAL

NOTICE! The air motor is **not** suitable for use with the Sirona DAC UNIVERSAL or with a cleaning and disinfection device.

We recommend using **Sirona DAC UNIVERSAL** for automated cleaning, disinfection, care, and unwrapped sterilization.

For further details, refer to the operating instructions supplied with the unit.

- ✓ The pushbutton chuck is manually maintained [→ 42].
 - ✓ The instrument is conditioned using the DAC UNIVERSAL.
1. Check whether the instrument is clean after conditioning under good lighting (min. 500 lux) and color rendering index (min. 80 Ra).
 2. If it is still dirty, repeat the process.

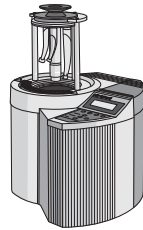



↪ For further conditioning, the instrument should be dry and free of residue.

3. Pack the instrument in packaging material suitable for sterilization and storage. e.g. paper/plastic composite packaging.
4. Sterilize the instrument [→ 43].

7.2.2 ... with cleaning and disinfection equipment

NOTICE! The air motor is **not** suitable for use with the Sirona DAC UNIVERSAL or with a cleaning and disinfection device.

The instrument can also be cleaned and disinfected in a suitable piece of cleaning and disinfection equipment.

The cleaning and disinfection equipment used must be approved by its manufacturer for the cleaning and disinfection of dental instruments and comply with EN ISO 15883-1 (e.g. 95°C (203°F) and 10 min. holding time).

7 Conditioning



36

Operating Instructions T4 LINE / Air Motor

For further details, refer to the operating instructions supplied with the unit.

- ✓ The instrument is conditioned with a cleaning and disinfection device.
- 1. Check whether the instrument is clean after conditioning under good lighting (min. 500 lux) and color rendering index (min. 80 Ra).
- 2. If it is still dirty, repeat the process.
 - ✎ For further conditioning, the instrument should be dry and free of residue.
- 3. Blow the instrument out with max. 3 bar.
- 4. Manual care of mechanical parts [→ 39].
- 5. Manual care of the push button chuck [→ 42].
- 6. Pack the instrument in packaging material suitable for sterilization and storage. e.g. paper/plastic composite packaging.
- 7. Sterilize the instrument [→ 43].

63 97 959 D3451.201.06.06.02

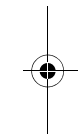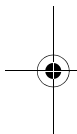

7.3 Manual cleaning and disinfection

IMPORTANT: Manual conditioning is possible in exceptional cases if the national/regional regulations are followed. The national/regional regulations are to be checked before.

NOTICE! Condition immediately, or at the latest, one hour after treatment.


NOTICE! Never clean in an ultrasound bath!

NOTICE! Never immerse in disinfectant solution!


NOTICE! Only use Sirona T1 Spray.

IMPORTANT: Use a soft, clean, and disinfected brush for cleaning.

- ✓ Wear appropriate protective clothing.
- ✓ All disinfectants must be approved in your country and have proven bactericidal, fungicidal and virucidal properties. Only use disinfectants with **no** protein-fixing properties.

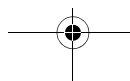


7 Conditioning

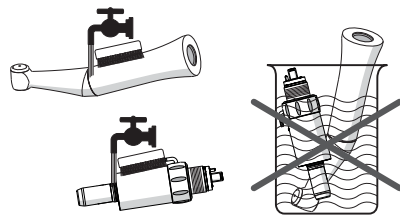

Operating Instructions T4 LINE / Air Motor

1. Brush T4 Line / Air Motor under running water ($< 38^{\circ}\text{C}$, $< 100^{\circ}\text{F}$, at least drinking water quality) and good lighting (min. 500 lux) and color rendering index (min. 80 Ra) until no more dirt can be seen, for at least 10 seconds.
2. Flush the drive channels with spray.
3. Conduct thermal disinfection or unwrapped steam sterilization.
4. Maintain mechanical parts manually [→ 39].
5. Maintain the push button chuck manually [→ 42].
6. Pack T4 Line / Air Motor in packaging material suitable for sterilization and storage. e.g. paper/plastic composite packaging.
7. Conduct sterilization [→ 43].

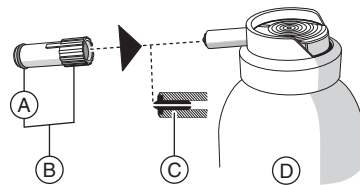

63 97 959 D3451.201.06.06.02

7.4 Manual maintenance

7.4.1 Maintenance of mechanical parts

Intervals


- At least every noon and evening
- Prior to each sterilization
- After every thermal disinfection without integrated maintenance

Required accessories

- | | |
|---|--------------------------|
| A | O-ring for spray adapter |
| B | Spray adapter |
| C | Black nozzle insert |
| D | T1 spray |

NOTICE! Only use Sirona T1 Spray.

7 Conditioning

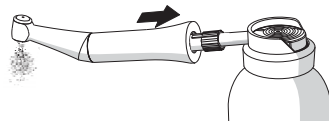


Process

Instrument

- ✓ The spray adapter is disinfected.
- ✓ The O-ring on the spray adapter is intact.
- 1. Fit the spray adapter onto the nozzle of the spray can.
- 2. Insert the instrument until it snaps into place and hold it.
- 3. Spray the instrument for 1 - 2 seconds.
IMPORTANT: Hold the spray can upright.
- 4. Wipe any spray that comes out with a disinfection cloth.
- 5. Repeat the process until the spray is clear.

Operating Instructions T4 LINE / Air Motor

Operating Instructions T4 LINE / Air Motor

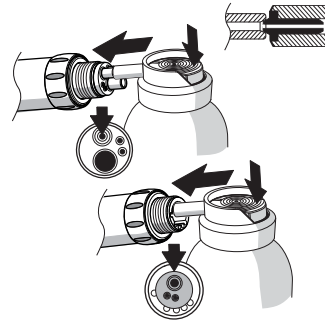

Air motor

✓ The instrument is removed [→ 23].

✓ The black nozzle insert is in the spray can tube.

1. Detach the motor from the supply hose [→ 22].

2. Set the drill direction to F.




3. Insert the nozzle of the spray can in the driving air tube. Hold the motor with a cloth.

4. Spray into the motor for 1-2 seconds.

IMPORTANT: Hold the spray can upright.

5. Wipe any spray that comes out with a disinfection cloth.

6. Repeat the process until the spray is clear.

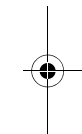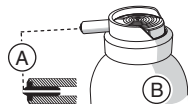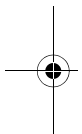

7.4.2 Care of the push button chuck

Use Sirona T1 Spray on the push button chuck to remove deposits and ensure proper functioning of the clamping system.

Interval



- At least once a week

Required accessories

A	Black nozzle insert
B	Sirona T1 Spray

NOTICE! Only use Sirona T1 Spray.

Operating Instructions T4 LINE / Air Motor

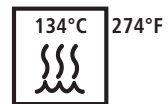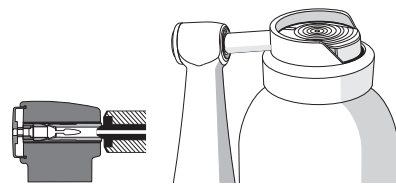

63 97 959 D3451.201.06.06.02

Procedure

- ✓ The instrument is clean and disinfected.
- 1. Press the contra-angle handpiece head with the chuck firmly against the spray can nozzle.
- 2. Spray the chuck for 1 - 2 seconds.
IMPORTANT: Hold the spray can upright.
- 3. Wipe any spray that comes out with a disinfection cloth.

7.5 Sterilizing

- ✓ T4 Line / Air Motor is cleaned and disinfected.
- ✓ T4 Line / Air Motor maintenance is complete.
- ✓ If necessary, T4 Line / Air Motor can be packed in packaging suitable for sterilization and storage, e.g. paper/plastic composite packaging or container.

7 Conditioning

7 Conditioning



44

Operating Instructions T4 LINE / Air Motor

- Sterilize T4 Line / Air Motor in the steam sterilizer with saturated water vapor.

Overpressure: 2.04 bar (29.59 psi)

Temperature: 134°C (274° F)

Holding time: 3 min.

Steam sterilizers that comply with either EN 13060 class B (for example, DAC PREMIUM / DAC PROFESSIONAL) or EN 13060 class S and are also suitable for the sterilization of straight/contra-angle handpieces are approved.

NOTICE! Do not exceed 140°C (284°F), even during the drying phase.

After sterilizing

1. Remove T4 Line / Air Motor from the steam sterilizer immediately.

⚠ CAUTION! T4 Line / Air Motor is hot.


63 97 959 D3451.201.06.06.02




6397959.book Seite 45 Mittwoch, 26. Juli 2017 2:52 14




Operating Instructions T4 LINE / Air Motor

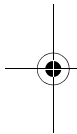

63 97 959 D3451.201.06.06.02

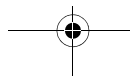

7 Conditioning

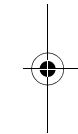

45



- NOTICE!** Do **not** try to accelerate the cool-down process by placing it in cold water. This will damage your T4 Line / Air Motor!
2. Store T4 Line / Air Motor so that it is protected against contamination.
 3. Sterilize again once the storage period has elapsed.

8 Maintenance

Elastomers, e.g. O-rings, must be replaced depending on their degree of wear.


Check the clamping system annually for safety reasons.


8.1 Replacing the sealing washer

If water leaks out between the motor and the hose connection, replace the sealing washer.

1. Disconnect the motor from the supply hose.
2. Use a probe etc. to remove the defective sealing washer (B) from the rear of the motor.
3. Attach the new sealing washer, observing the position of the tubes and contact pins and press it to the stop.

IMPORTANT: The bead (A) on the sealing washer faces the supply hose.

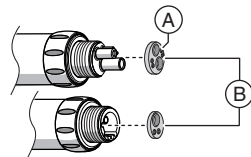

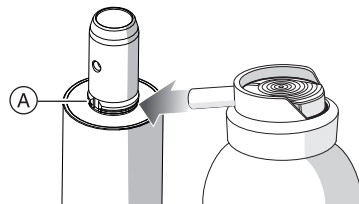

8.2 Lubricate the lock washer

Interval

- 1x weekly

Procedure

1. Spray some T1 spray on the stop spring ring (A).
2. Turn the stop spring ring to distribute the T1 spray.

9 Spare parts and consumables

Use only original Sirona parts or parts approved by Sirona.

	REF		REF
Pin for contra-angle handpiece bur	89 17 866	Tool for pin (contra-angle handpiece bur)	89 17 874
NaCl silicone hose, thin, 200 mm long, 6 pcs.	59 17 419	Spray clip B for contra-angle handpiece	64 28 374
Midwest gasket	58 55 460	Borden gasket	59 50 121
O-ring 8 x 1	70 36 189	T1 spray (6 x 250 ml cans)	59 01 665
Spray adapter	89 17 858	O-ring for spray adapter	70 36 353


10 Storage and transport conditions




Protect from
moisture




Sensitive
contents

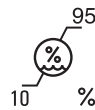

Relative
humidity


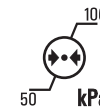


Air pressure

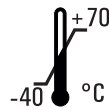


Temperature








After a severe change in
temperature, allow sufficient
time for acclimation.

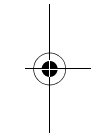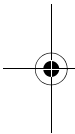

11 Disposal

- According to current information, the product does not contain any substances that are hazardous to the environment.
- Disinfect the product prior to disposal.
- Observe the applicable disposal regulations for your area.

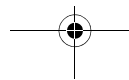

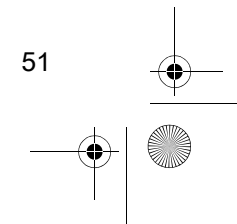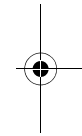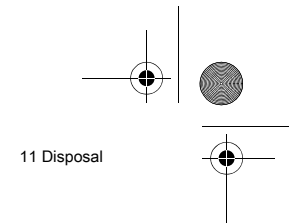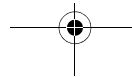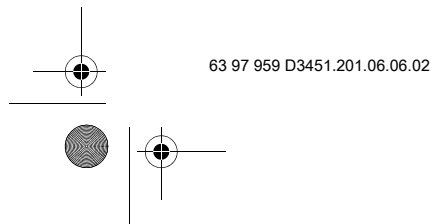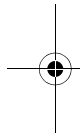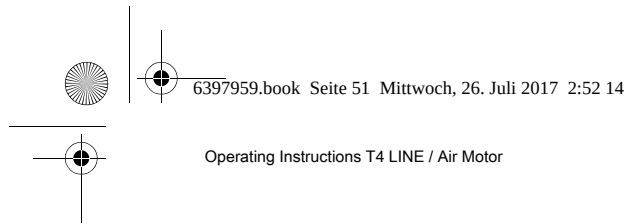

We reserve the right to make any alterations which may be required due to technical improvements.

D3451.201.06.06.02 07.2017

Ä.-Nr.: 124 093

CE
0123

 **Sirona Dental Systems GmbH**
Fabrikstraße 31
D-64625 Bensheim
Germany
www.sirona.com

in USA: Sirona Dental Systems, LLC
4835 Sirona Drive, Suite 100
Charlotte, NC 28273
USA
www.sirona.com

63 97 959 D3451